


Pressure Area Care
POLYFLOAT™ Dormira Mattresses
...safe, effective, harm free care


Pressure Reducing Reactive Therapy

The POLYFLOAT Dormira mattress is designed to *redistribute interface pressure, reduce shear and provide comfort* on profiling or non-profiling beds.

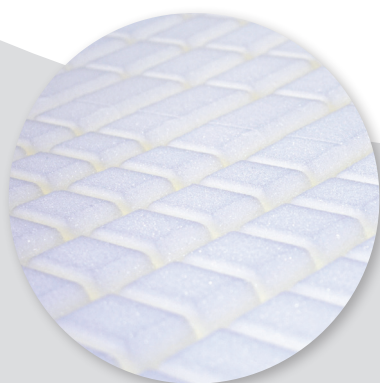


OPTIMISED IMMERSION AND ENVELOPMENT

The POLYFLOAT Dormira features a series of castellations across the mattress surface. These assist the foam to conform to individual body contours, which redistributes pressure by allowing the partial immersion and envelopment of the patient.

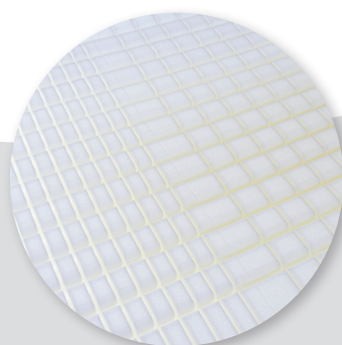

ZONED MATTRESS

The POLYFLOAT Dormira mattress has dedicated head, torso and heel zones to provide the patient with optimal support and pressure redistribution across the body.

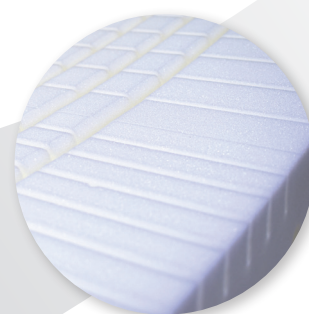


ASSISTED PATIENT TRANSFER

Intentionally wider edges on each side of the POLYFLOAT Dormira mattress enhance patient stability and assist with patient transfers and care.

Polyfloat™

Dormira

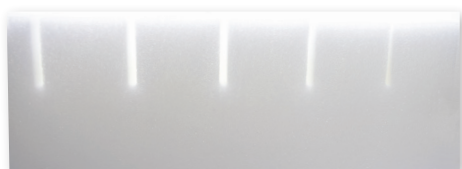

TECHNICAL DATA

Length:	1980 mm
Width:	880 mm
Depth:	150 mm
Weight:	11 kg
Weight limit:	185 kg
Foam density:	34kg/m ³
Fatigue rating:	Very severe



MULTI-STRETCH COVER

The multi-stretch cover reduces shear and friction during patient movement or manual repositioning. The material is waterproof and moisture-vapour permeable, contributing to the management of the skin 'microclimate' between the patient and support surface. All seams are welded to protect the inside of the mattress from any fluid ingress.

POLYFLOAT Dormira is a Trade Mark of Talley Group Limited.

Talley manufacture products to comply with National and International safety standards and are certified to ISO 13485 and Directive 93/42/EEC. Every care has been taken to ensure that the information contained in this brochure was correct at the time of going to press. However, Talley reserve the right to modify the specification of any product without prior notice in line with a policy of continual product development. Our standard terms and conditions apply. © Talley Group Limited 2007 - 2018. All rights reserved.



Talley Group Limited
Premier Way, Abbey Park Industrial Estate
Romsey, Hampshire, SO51 9DQ England
TEL: +44(0)1794 503500 FAX: +44(0)1794 503555
EMAIL: sales@talleygroup.com

www.talleygroup.com

